BÁSCULA DIGITAL DE CRISTAL

Manual de instrucciones

MED-095

V0.0 0717m

IMPORTANTE

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V0.0 0717m

IMPORTANT


Read these manual to avoid any damage.
The instructions of these manual are for reference about the product.
There may be differences due to update.
Please check our web site www.steren.com to obtain the latest version of the instruction manual.

CAUTIONS

- Always use a soft and dry cloth to clean the scale.
- Do not use the scale if it is wet or if you have wet feet.
- If you will not be using it for an extended period of time, remove the batteries.
- Always place the scale on a level surface.
- Do not exceed the maximum capacity (396 lb).
- Avoid device dropping, as it may be damaged.
- This product cannot be used by people with different abilities, unless they have preparation and supervision.
- This product is NOT a toy; keep it away from children.

BEFORE TO USE

- 1 Remove the cover located in the rear part of the scale. Place the **2 AAA batteries** correctly. *Immediately the screen will turn on*

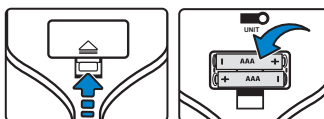- 2 Select the weight unit (**kg, lb or st:lb**) pressing the button as shown in the below image.

HOW TO USE

- 1 Place both feet on the scale and keep as still as possible. *The screen will turn on immediately.*
- 2 Wait for a few seconds while the reading is stabilized.
- 3 *The final reading will appear flashing on the screen. Wait for a few seconds until it is fixed. At finish you can get off from the scale.*

SPECIFICATIONS

Input: 3 V --- (2 x AAA)